# **Infrared Spectroscopy**

K. Sarath D. Perera


Senior Professor in Chemistry, Department of Chemistry, The Open University of Sri Lanka

The energy of most of the molecular vibrations falls into the infrared region (4000 - 400 cm<sup>-1</sup>) of the electromagnetic spectrum. Infrared (IR) or vibrational spectroscopy is used by organic chemists to: (i) identify functional groups in solid, liquid, or gaseous forms, (ii) differentiate inter- and intra-molecular hydrogenbonding, and (iii) investigate chemical reactions. Recent developments in Fourier Transform Infrared (FT-IR) spectroscopy extended its applications to atmospheric, archaeological, pharmaceutical, food and biological fields, and also to analyze compounds/materials qualitatively and quantitatively.

IR spectrometer measures the interaction of infrared radiation with matter by means of absorption/ transmittance. When molecules absorb energy, they get excited to higher vibrational/rotational energy levels, which causes constant movements in each covalent bond with respect to other atoms present in the molecule. These movements are called fundamental vibrational modes (FVM).


#### Fundamental vibrational modes

There are two main FVM: stretching and bending. **Stretching** modes are of two types: *e.g.*, symmetric and asymmetric.

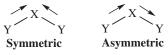


**Bending** modes are of two types: *e.g.*, (i) in- plane (rocking and scissoring) and (ii) out of plane (twisting and wagging).

#### In plane bending

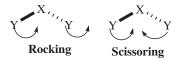

# Out of plane bending

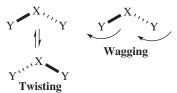

The number of **vibrational modes** (degrees of freedom) for a **linear** molecule is equal to 3N-5 while that for **non-linear** molecule is 3N-6, where N = number of atoms in the molecule. For example, the number of modes in HCl and  $CO_2$  are 1 and 4, respectively. The non-linear molecule, benzene has 30 modes.

# Infrared spectra

An infrared spectrum is a graph of percentage transmittance (%T) on the vertical axis versus wavenumber ( $\tilde{v}$ ) in cm<sup>-1</sup> on the horizontal axis (Fig. 1).

$$%T = \frac{\text{transmitted intensity x 1}00}{\text{original intensity}}$$

If the molecule does not absorb radiation, then %T is 100% and no peaks appear in the spectrum. When a molecule absorbs energy the %T decreases. The wavenumber range is related to the bond type and the **functional group (FG)** present in the molecule. Not all covalent bonds display IR bands; only the polar bonds are **IR active**. The **intensity** of the band depends on the magnitude of the **dipolar moment**. For example, polar bonds such as carbonyl (C=O) groups produce **strong** (s) bands; asymmetric stretching of bonds display **medium** (m) bands while symmetric stretching of bonds show **weak** (w) bands. Overlapping of a weak band on a side of a strong band gives rise to a **shoulder** (sh). The vibrational frequency of a bond (v or v) can be calculated using Hooke's law,

$$v = \frac{1}{2\pi} \sqrt{\frac{K}{\mu}}$$
 or  $\bar{v} = \frac{1}{2\pi c} \sqrt{\frac{K}{\mu}}$ 

K = force constant, c = velocity of light, reduced mass =  $\mu = (m_1 m_2)/(m_1 + m_2)$ ,  $m_1$  and  $m_2$  are masses of atoms, strong bonds show higher frequencies,  $\nu = c\bar{\nu}$ , c = speed of light.

The IR spectrum is divided into two parts:

- Functional group region (4000-1400 cm<sup>-1</sup>)
- Finger print region (1400-400 cm<sup>-1</sup>) which is specific to a molecule (see Fig. 1)

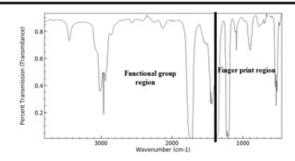


Figure1: Two regions of an IR spectrum

# Characteristic IR frequencies of FGs

Even a simple organic molecule can have a large number of vibrational modes, *e.g.*, C-H, C-C, C-O, O-H, N-H, C=O, C=N *etc.* Bending and stretching frequencies of organic molecules are given in Table 1. The following bands are suggested by comparing the values obtained for a large number of molecules. RNH $_2$  and RCONH $_2$  show two strong bands at ~3340 & 3300; R $_2$ NH gives only one band. The v(C=O) values (1870-1650 cm $^{-1}$ , see Table 1) depend on the inductive and resonance effects exerted by the atoms/groups attached to the C=O group.

**Table 1:** Bending & stretching frequencies of common organic molecules in cm<sup>-1</sup>


Type of molecules	Bond	Bending frequency	Stretching frequency
Alkane	С-Н	1470-1450	3000-2840m
Alkene	C=C		1680-1620
	=С-Н	1000-650s	3090-3010m
Alkyne	C≡C		2260-2100w
	≡С-Н	700-600	3330-3270
Aromatic	C=C		1600-1585 1500-1400
	=С-Н	900-690s	3100-3010m
Anhydride			1870-1800 I 1820-1770 II
RCOCl			1820-1770
Ester	C=O		1820-1770
RCO <sub>2</sub> H			1725-1700
Aldehyde			1750-1720s
Ketone			1730-1705s
Amide	]		1690-1620

Alcohol	C-O		1260-1050
			3500-3200
H-bonded	О-Н	1620-1590	3600-3200
Carboxylic acid		1440-1395 950-940	3400-2400
Amine 1 <sup>ry</sup>	N-H	1650-1580	3400-3300w 3330-3250w
Amine 2ry	N-H		3350-3310w
	C-N		1360-1180s
Amide	N-H	1650-1560	3370 (asym) 3150 (sym)
Nitrile	C≡N		2280-2210m
Imine	C=N		1690-1640m
R-NO <sub>2</sub>	N=O		1570-1500s 1385-1300s

Note that v(asymmetric) > v(symmetric).

Absorption bands for the C-X group of alkyl halides appear in the finger print region, e.g.,  $\nu$ (C-F)  $\sim$  1200 cm<sup>-1</sup>,  $\nu$ (C-Cl) 850-550 cm<sup>-1</sup>,  $\nu$ (C-Br) 690-510 and  $\nu$ (C-I) < 650 cm<sup>-1</sup>.

The summary of absorption bands due to molecular vibrations is given is Fig, 2.

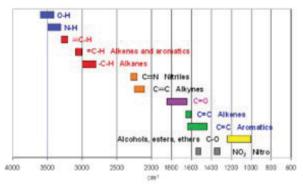

**Figure 2:** Summary of absorption bands due to molecular vibrations

# Factors affecting the IR frequency

The main factors that affect the IR frequencies of a molecule are conjugation effect, resonance effect, inductive effect {electron donating (+I) and withdrawing (-I) effects}, hydrogen bonding, the ring size of the molecule, *etc*.

Vol. 39 No. 1, June 2022

#### a) Conjugation and resonance effects

The resonance effect of the vinyl group attached to methylvinyl ketone increases the single bond character of the C=O group, thus, the  $\nu$ (C=O) value is lowered to 1680 cm<sup>-1</sup> (Fig. 3).

**Figure 3:** Resonance effect on the C=O group upon extended conjugation

Similarly,  $\nu(C=O)$  shifts to lower values when an aryl group is conjugated to it, *e.g.*, the  $\nu(C=O)$  values for PhCHO and PhCO<sub>2</sub>H are 1700 and 1680 cm<sup>-1</sup>, respectively.

#### b) Inductive effect

The neighboring groups attached to the particular functional group strongly affect the IR frequency of that functional group. The trends in  $\nu(C=O)$  for some carbonyl compounds are given in Fig. 4.

**Figure 4:**  $\nu$ (C=O) values of C=O groups

In chloro-derivatives of acetone, the electron withdrawing effect of chlorine increases, thus, the bond order of the C=O group,  $\nu$ (C=O), increases. The bond order of the C=O group in formaldehyde is higher than that of acetone. The donation of the lone-pair electrons on amino nitrogen weakens the C=O bond. The amino group can have a combination of inductive and resonance effects. Electron withdrawing effect of electronegative atoms such as oxygen and fluorine increases the bond order of the C=O group.


# c) Hydrogen bonding

68

Inter- and intra-molecular H bonding (Fig. 5) lowers the frequency of the O-H bond. Neat phenol gives a

broad band around 3350 cm $^{-1}$  for the O-H bond which is H bonded. Upon dilution, the  $\nu$ (O-H) value shifts towards a higher value,  $\sim$ 3600 cm $^{-1}$ .

inter-molecular H bonding intra-molecular H bonding

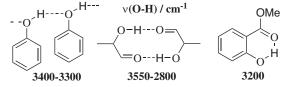

**Figure 5:** Inter- & intra-molecular H-bonding of O-H bonds

The IR spectrum of acetic acid (Fig. 6) showed a broad band between 3550-2800 cm<sup>-1</sup> for the hydrogen bonded OH group.

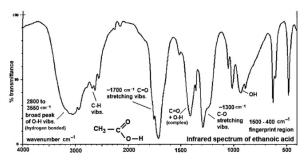


Figure 6: IR spectrum of acetic acid

Note that stretching frequencies  $\nu$ (C=O) and  $\nu$ (C-O) of the carboxylic group are around 1700 and 1300 cm<sup>-1</sup>, respectively.

Hydrogen bonds are not present in the gas phase. The IR spectrum of the **gaseous** HCHO (Fig. 7) showed sharp bands at 2850 (asymmetric), 2750 (symmetric), 1485 (scissoring), 1250 (rocking) and 1165 (wagging) as the vibrational modes of the  $CH_2$  group. The stretching frequency of  $\nu(C=O)$  is 1750 cm<sup>-1</sup>.

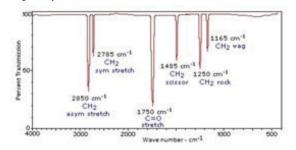



Figure 7: IR spectrum of HCHO vapor

## d) Ring size effect

The  $\nu(C=O)$  values given in Table 1 are normally for unstrained acyclic carbonyl compounds. The

Vol. 39 No. 1, June 2022

vibrational frequency of a C=O group in a strained ring system increases as shown in Fig. 8. Acyclic acetone and cyclohexanone have the same value of 1715 cm<sup>-1</sup>. The  $\nu$ (C=O) value gradually increases as the size of the ring decreases due to increase in the angle strain. Similar trend is observed for strained cyclic esters or lactones.

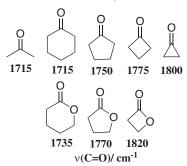

**Figure 8:** Ring size effect in cyclic ketones & lactones

# **Metal Carbonyls**

IR spectroscopy is an important tool to study metal carbonyls because the C $\equiv$ O stretching frequencies  $\nu$ (C $\equiv$ O) give rise to **strong** bands between 2150 and 1600 cm $^{-1}$ . The free CO has a  $\nu$ (C $\equiv$ O) value of 2143 cm $^{-1}$ .

IR spectra of metal carbonyls provide the following information:

- Bonding mode of CO whether it is terminal, doubly bridging or triply bridging.
  - The  $v(C\equiv O)$  for terminal and bridging carbonyls decreases in the following order:
  - M-CO (2150-1900) >  $M_2(\mu_2$ -CO) (1900-1750) >  $M_3(\mu_3$ -CO) (1750-1600) in cm $^{-1}$ .
- Geometry of the molecule the number of bands indicates the symmetry of the M(CO)<sub>n</sub> fragment.
  Metal carbonyls, [Ni(CO)<sub>4</sub>], [M(CO)<sub>6</sub>], trans-[M(CO)<sub>4</sub>L<sub>2</sub>] (M = Cr, Mo, W), and [M(CO)<sub>3</sub>(L<sub>axial</sub>)<sub>2</sub>] (M = Fe, Ru) show only **one** IR band whilst fac-[M(CO)<sub>4</sub>L<sub>3</sub>] (M = W, Mo) shows **two** IR bands.
- 3. Electron donor and withdrawing nature of other co-ligands electron donors lower the stretching frequency of C≡O.

Complex	$\nu(CO)/cm^{-1}$
$fac$ - $[Mo(CO)_3(PCl_3)_3]$	2040, 1991
$fac$ - $[Mo(CO)_3(PEt_3)_3]$	1937, 1841

4. Oxidation state of the metal or the charge on the complex ion - anionic complexes show lower stretching frequencies.

Complex/ion	$v(CO)/cm^{-1}$
$[\mathrm{Mn(CO)}_6]^{\scriptscriptstyle +}$	2090
$[Cr(CO)_6]$	2000
$[V(CO)_6]^-$	1860

### Metal hydrides

The IR absorption frequency of the M-H bond,  $\nu$ (M-H), appears in the region of 1500-2200 cm<sup>-1</sup>, but intensities are often **weak** to **medium** when compared to intensities of the carbonyl bands of metal carbonyls.

The presence of a M-H bond can be confirmed by converting it into a M-D bond with  $D_2O$  and recording the IR spectrum again. The absence of the M-H band and appearing of a new M-D band is observed according to the relationship,  $v(M-D) = v(M H) \div \sqrt{2}$ .

### Other metal complexes

It is important to have an idea of the values of stretching frequencies of metal complexes and neutral ligands, used in coordination and organometallic chemistry. Table 2 gives the ranges of stretching frequencies of some neutral ligands (L) when they are **bound** to a metal center or **unbound**, *i.e.*, the **free** ligand.

**Table 2:** Stretching frequencies of some neutral ligands (L) in cm<sup>-1</sup>.

L	ν(free-L)	ν( <b>bound-</b> L)
C≡S	1270	1410-1160 (terminal)
		1150-1020 (μ-CS)
N≡N	2330	2200-1800 (terminal)
		2030-1700 M <sub>2</sub> (μ-N <sub>2</sub> )
N≡O	1880	1830-1610 (linear)
		1720-1520 (bent-NO)
N≡O+	2250	1520-1490 M <sub>2</sub> (μ-NO)
CH <sub>2</sub> =CH <sub>2</sub>	1623	1520-1450 M(η²-L)
PhC≡CH	2110	2200-1900 (MC≡CPh)
		2200-1900 (linear)
MeN≡C	2170	1870-1830 (bent)
		1700-1600 M <sub>2</sub> (μ-L)
Н-Н		2900-2300 M-(η <sup>2</sup> -H <sub>2</sub> )

# Inorganic ions

Table 3: Stretching frequencies of common inorganic ions in  ${\rm cm}^{-1}$ .

Ion	ν(ion)/cm <sup>-1</sup>
NH <sub>4</sub> <sup>+</sup>	3300-3030, 1485-1390
CO <sub>3</sub> <sup>2-</sup>	1450-1410, 880-800
SO <sub>4</sub> <sup>2-</sup>	1130-1080, 680-610
NO <sub>3</sub> -	1380-1350, 860-800
NO <sub>2</sub> -	1250-1230
PO <sub>4</sub> <sup>3-</sup>	1100-1000
NC <sup>-</sup> , NCS <sup>-</sup> , SCN <sup>-</sup>	2200-2000
SiO <sub>4</sub> <sup>2-</sup>	1100-900
MnO <sub>4</sub>	920-890, 850-840

- -